



A NEW TWIST ON PILL CRUSHING



MAXGRIND[™] Pill Crusher & Grinder

Eloquently designed with a dynamic crushing and grinding action to easily deliver finely crushed medication from the top down.

www.maxpertmedical.com • 866.231.1222
info@maxpertmedical.com

CRUSHER. GRINDER. PILL STORAGE. SERVING CUP.




MAXGRIND™ is a revolutionary, high-quality twist pill crusher and grinder featuring an ergonomic, easy-to-grip and twist design.

There is no need to unscrew the top to check progress or retrieve medications! Retrieving crushed medication is simplified by the top down process where crushed medication is sifted and funneled through patterned slots into the clear, easy-view bottom serving cup.

- **Comfortable anti-slip rubber grip.**
- **Easy dual action crushing and grinding.**
- **Exclusive top down crushing into serving cup.**
- **Free spin, pepper grinder action, delivers fine powder.**
- **Safe! BPA, latex-free and FDA compliant materials.**
- **Built-in storage with hinged lid.**

THE FIRST OF ITS KIND hand-held twist pill crusher that allows medication to be delivered in two ways:

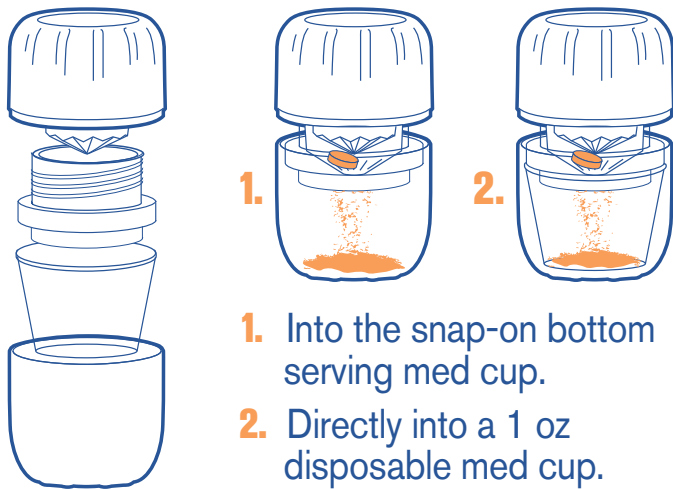

1. Into the snap-on bottom serving med cup.
2. Directly into a 1 oz disposable med cup.



Latex Free

Ideal for the every day home user or the professional providing patient specific pill crushing.

**CRUSHING
FROM
THE
TOP
DOWN**

Crush into the serving cup and either empty or mix and use as a mini bowl or directly into the disposable med cup to empty or mix and consume with little cleanup remaining.

1 step!

EASY DOES IT! FREE SPIN WITH LITTLE EFFORT.



FREE SPIN – Unlike traditional twist crushers which stop crushing when the top is screwed tight against the bottom, MAXGRIND features a unique screw design which continues to grind medication when tightened. Once the top is screwed down to initially break and crush pills, a simple reverse twist then slight pressure on the top allows the grinding surface to freely spin similar to a pepper grinder, continuing to grind medications without any resistance and with total ease.



Internal Cleaning

Feature - If you continue to free spin in the upright position after crushing is complete, air circulates within the crusher blowing remaining powder into the bottom cup, resulting in minimal medication remains on the grinding surface and in the crushing well.



Clear, snap on serving cup allows easy viewing and collection of crushed medication. Cup can be used as a serving cup to mix in and consume from directly or the crushed medication can be easily emptied over into food or drink or another container. Easy to mix and use with a normal spoon.

Ridges on the grinding surface provide fast, dual action crushing and grinding with minimal residue left behind.



Uniquely shaped and spaced ridges

on the top crushing surface both crush and grind medications into a fine powder that is sifted/funneled into the bottom serving cup. These ridges are designed so that powder is pushed out and funneled through the crusher well slots into the bottom serving cup, leaving less residue behind which ensures a more complete dose.

MAXGRIND™ also provides a simple change of tone when crushing is completed, serving as an indicator that crushing is completed.



PO Box 111
Galloway, OH 43119



MAXGRIND™ Pill Crusher & Grinder

Top Pill Storage Compartment with Hinged Lid - Great for storing pills when traveling.

Comfortable and easy to grip anti-slip symmetrical design - Ergonomic design and easy to grip rubberized material provide comfort when crushing and aid in the twisting effort.

Easy View Snap-On Bottom Serving Cup - Can be used as a mini bowl.



Multiple Pills Can Be Crushed At One Time:

MAXGRIND's smart crushing surface and design allows crushing and grinding of multiple pills at once with minimal effort.

Easy Maintenance: After the internal cleaning feature helps release maximum powder into the medication cup, **MAXGRIND** can be rinsed with water for easy cleaning and preparation for next use.

Dimensions: 3-3/8"L x 2-3/16"W

Sold individually: Available in Red or Blue.

Reorder: Blue: BMG01B Red: BMG01R



MAXGRIND™ is a trademark of Maxpert Medical, LLC.

*Patent pending